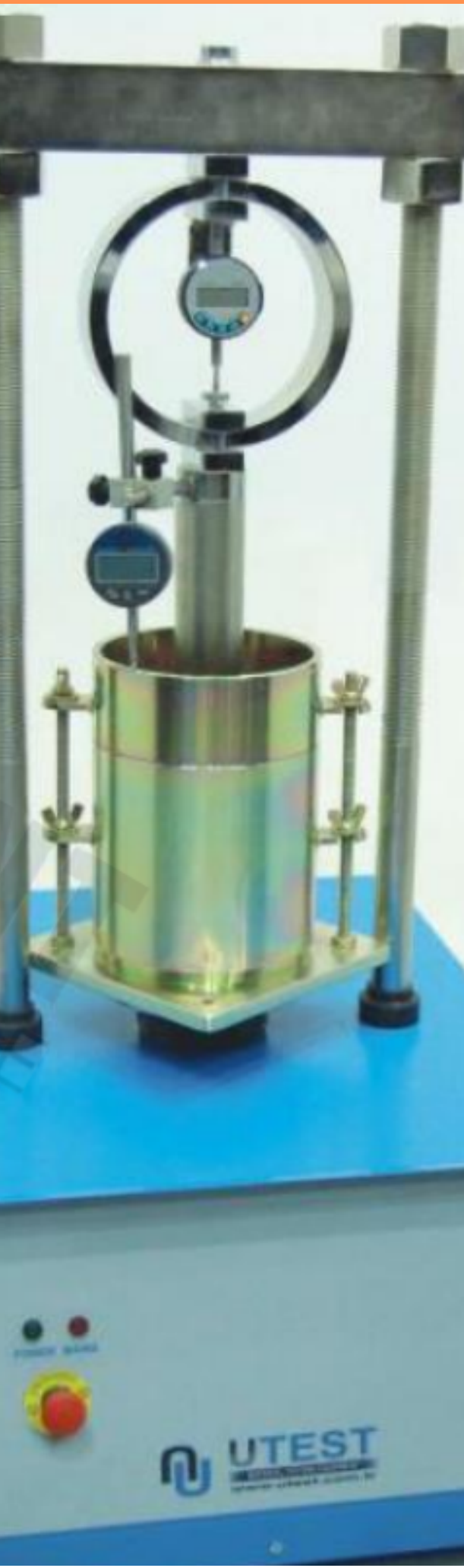




• Product Code

UTS-0852 CBR Test Machine with Load Ring, 50kN capacity

Models for 220-240V 50-60 Hz, 1 ph.

UTS-0852

Models for 110-120V 60 Hz, 1 ph.

UTS-0852-N

• Standards

EN 13286-47; BS 1377:4; ASTM D1883; AASHTO T193; NF P94-078; AS 1289.6.1.1; UNI CNR 10009

The UTS-0852 is designed for performing laboratory evaluation of the CBR value of highway sub bases and subgrade and for the determination of strength of cohesive materials which have maximum particle sizes less than 19mm(3/4").

The UTS-0852 is designed to load the penetration piston into the soil sample at a constant rate to measure the applied load and piston penetration at predetermined intervals.

The machine is designed to be mounted on a suitable bench and comprises of a robust and compact two column frame with adjustable upper cross beam. The frame has 50kN capacity. The test speed is 1.27mm/min. for ASTM/EN/AASHTO/BS/NF tests. Loading and unloading are down from the front panel by UP/DOWN buttons. Unloading speed is adjusted 5mm/min for easy re-testing.

The CBR Test Machine is supplied complete with;

- Load Ring, 50 kN
- Digital dial gauge 25 x 0.01 mm (UTGM-0148) with holder (UTS-0853)
- Penetration Piston (UTS-0870)

Dimensions

480x650x1150 mm

Weight (approx.)

110 kg

Power

370 W